


# Safe Mitten II

**Mitten-style gloves designed to limit grasping for self-harm prevention.**

**Features "Idea Hooks" that are difficult for the wearer to open.**

**Breathable side mesh and cool-touch inner lining to reduce stuffiness.**

**Front zipper for quick access to hands for care or inspection.**

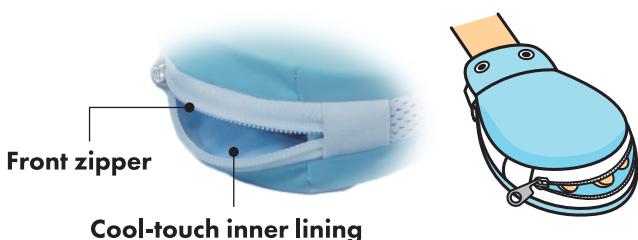

**Allows visual inspection of the inside of the gloves.**

**Allows hand inspection without removing the gloves.**

**Features a zipper with a locking function.**



※Please use the included special zipper pulls for opening and closing.



This product is used temporarily to support medical and caregiving procedures when there is a risk of the patient engaging in dangerous behaviors, such as self-harm or pulling out catheters, or when they try to tamper with their diaper.

## "Idea Hook"



**[How to unhook]**  
Support the back of the hook with your finger, and press the center with your fingernail.



**[How to close]**  
Place both hands on the edges of the hook over the flap and press firmly until you hear a 'click'.

## Material

Glove fabric : Nylon  
Inner fabric : Polyester  
Side mesh : Polyester  
Hoop-and-loop fastener : Nylon  
"Idea Hooks" : Stainless steel  
Wrist Strap : Polyester  
Protective Pad : Polyethylene

## Washing Instructions

Make sure hooks and buckles are closed before washing.  
Washing them open or detached can lead to wear or damage.



## Size

Wrist	13~22cm
Hand Length	14~20cm



## Specification

Product Code	Item	Unit
041051	Safe Mitten II	1 Pair
041053	Replacement Protective Pad	1 Piece
041055	Replacement Flap with "Idea Hooks"	1 Piece



**Taketora®**

<https://www.taketora-web.com/>

**Taketora Co., Ltd.**

[Human Care Department]

8F Yokohama Media Business Center, 2-23 Otamachi,

Naka-ku, Yokohama-shi, Kanagawa 231-0011, Japan

TEL +81-45-680-0390 FAX +81-45-680-0399